


Separation of 5 Agricultural Chemicals – 5

Data# E-034

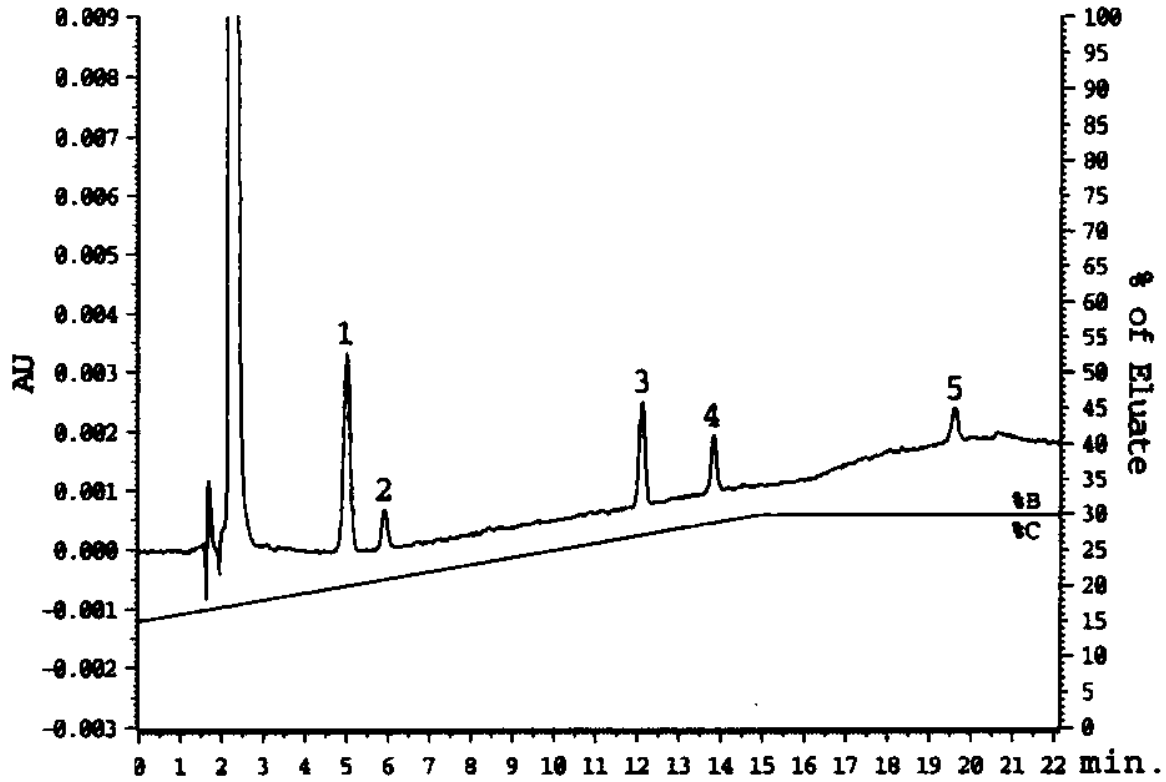

Column	Mightysil RP-18 GP 150-4.6mmØ 5µm	
Eluent	A; Water, B; Acetonitrile C; Methanol	
Flow Rate	1.0mL / min.	
Detection	UV270nm	
Temperature	40°C	
Sample Volume	40µl	
Detected peaks	1. Imidacloprid	0.1 µg / ml
	2. Benomyl	0.1 µg / ml
	3. Thiophanate-methyl	0.1 µg / ml
	4. Carbaryl	0.1 µg / ml
	5. Linuron	0.1 µg / ml