


Mightysil Application DATA

Separation of Pesticides standard solution 28

Data# E-040

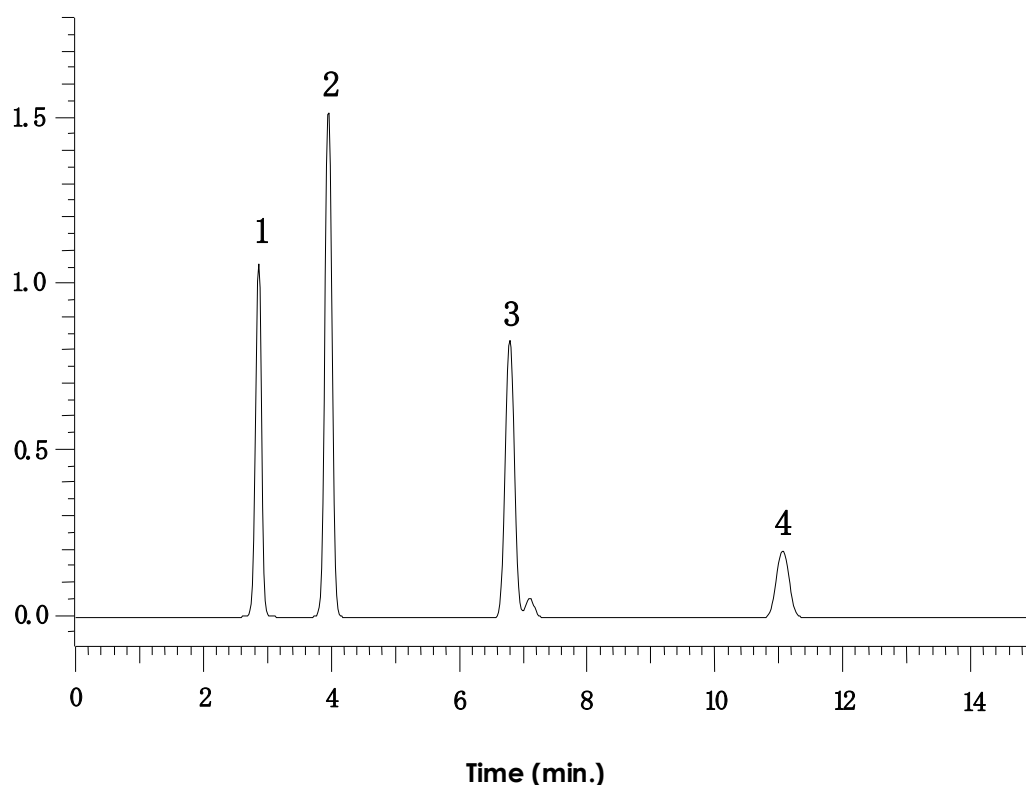


Column	Mightysil RP-18PA 250-4.6mmØ 5µm	
Eluent	Acetonitrile/ 50mM KH ₂ PO ₄ = 55/ 45 PH adjust: addition of Phosphoric acid to pH3.3	
Flow Rate	1.0mL/min.	
Detection	UV270nm Asulam UV230nm thiophanate-methyl, siduron, iprodion	
Temperature	40°C	
Sample Vol.	20µL	
Samples	1. Asulam	50ppm
	2. Thiophanate-methyl	100ppm
	3. Siduron	100ppm
	4. Iprodion	50ppm




Kanto Kagaku

Mightysil Application DATA


1

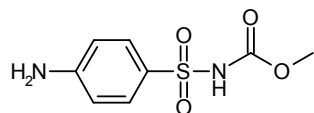

2

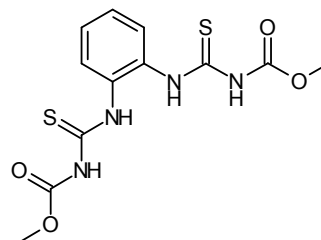

3

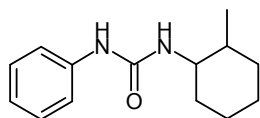

4

