


## Separation of Drugs for Animals – 2

Data# F-015

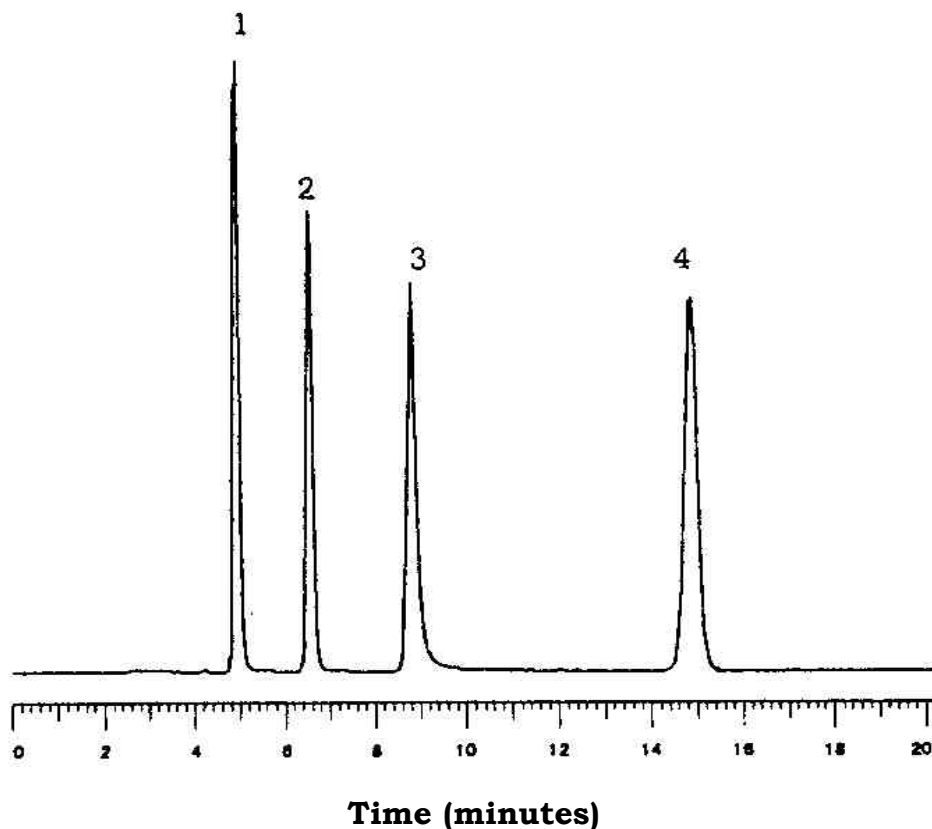

Column	Mightysil RP-18 GP 150-4.6mmØ 5µm
Eluent	Acetonitrile/Phosphate buffer pH3.0 = 15/85
Flow Rate	0.5mL / min.
Detection	UV285nm 0.03AUFS
Temperature	40°C
Sample Volume	5µL
Samples	1. 5-Hydroxythiabendazole 5µg/mL 2. 5-propylsulfonyl-1H-benzimidazole-2-amine 5µg/mL 3. 2-Quinoxalinecarboxylic acid 15µg/mL 4. Sulfadimidine 7µg/mL