


Comparison of Performance of Reversed Phase Packing Materials-3

Data# S-015

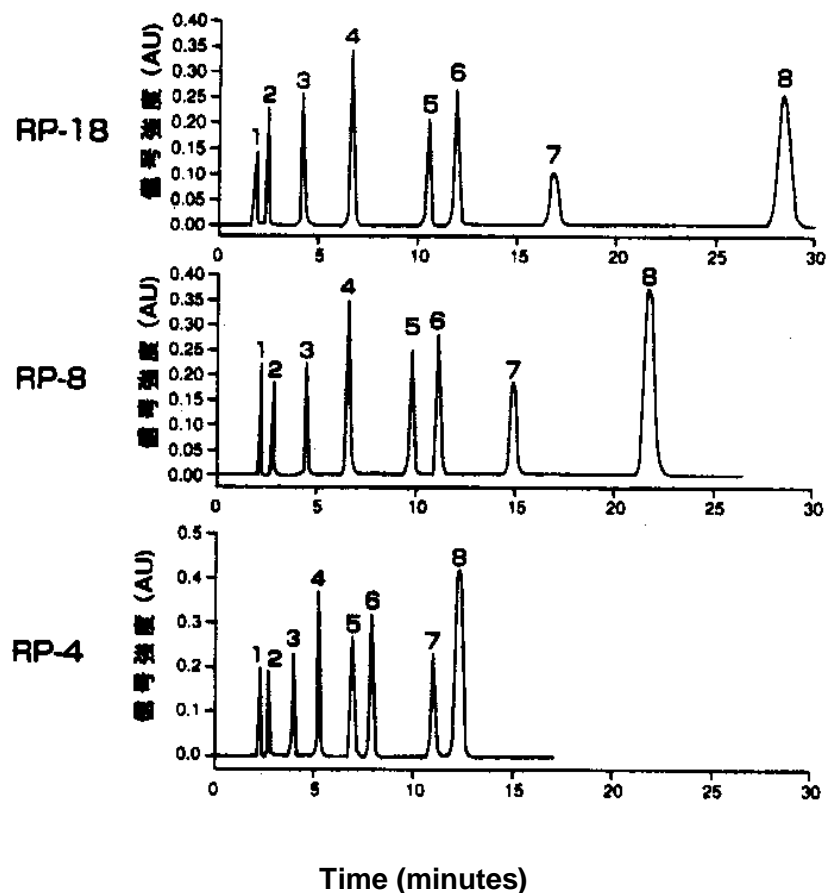

Column	Mightysil GP 150-4.6mmØ 5µm	
Eluent	Acetonitrile/Water = 30/70	
Flow Rate	1.0mL / min.	
Detection	UV254nm	
Temperature	40°C	
Sample Vol.	10µL	
Samples	1. Caffeine 35µg/mL 2. Pyridine 35 µg/mL 3. Aniline 85 µg/mL 4. o-Toluidine 250 µg/mL	5. N-Methylaniline 50 µg/mL 6. 2-Ethylaniline 310 µg/mL 7. 3-Nitroanisoole 140 µg/mL 8. N,N-Dimethylaniline 140 µg/mL